




Mastectomy – C-Frame





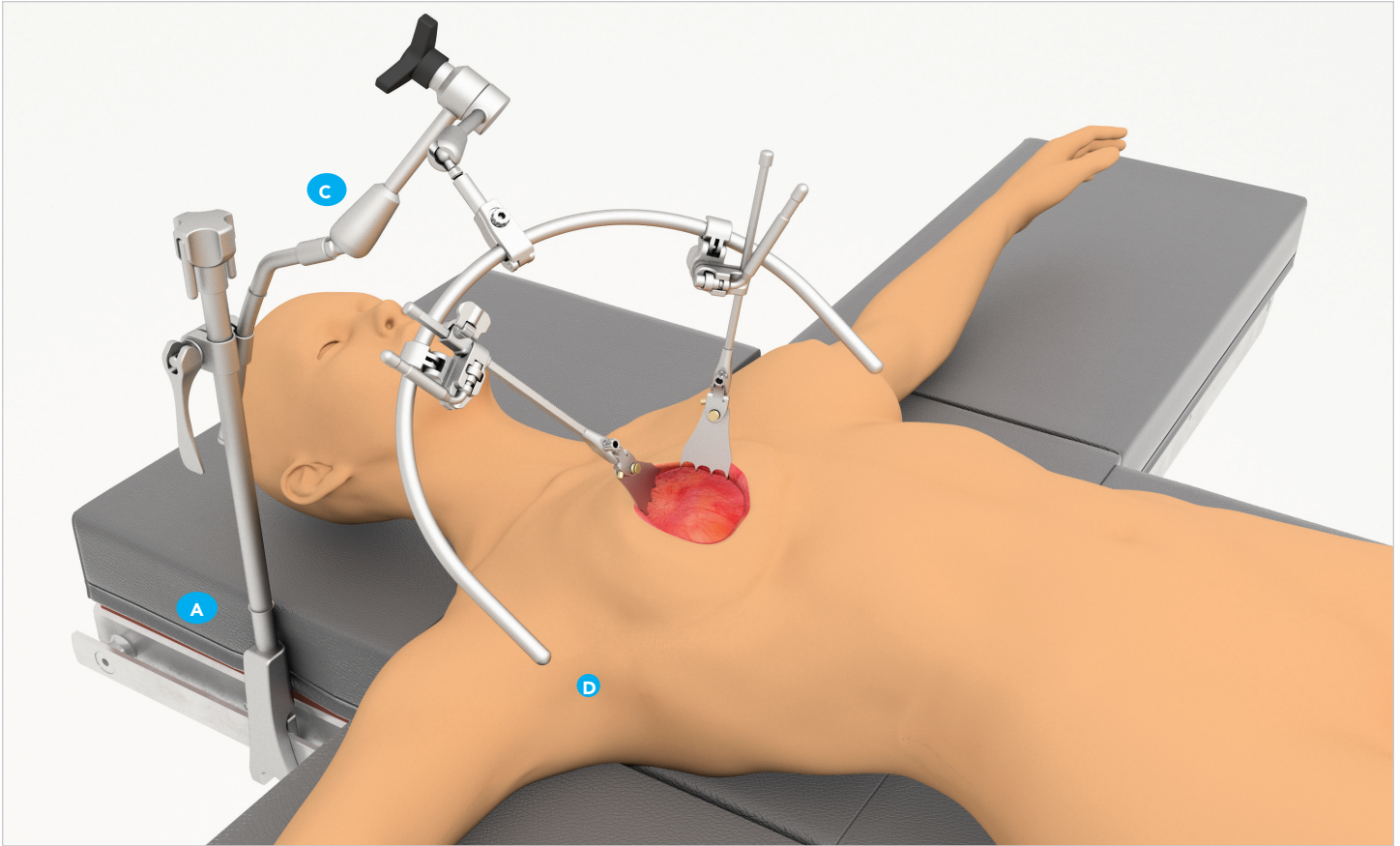

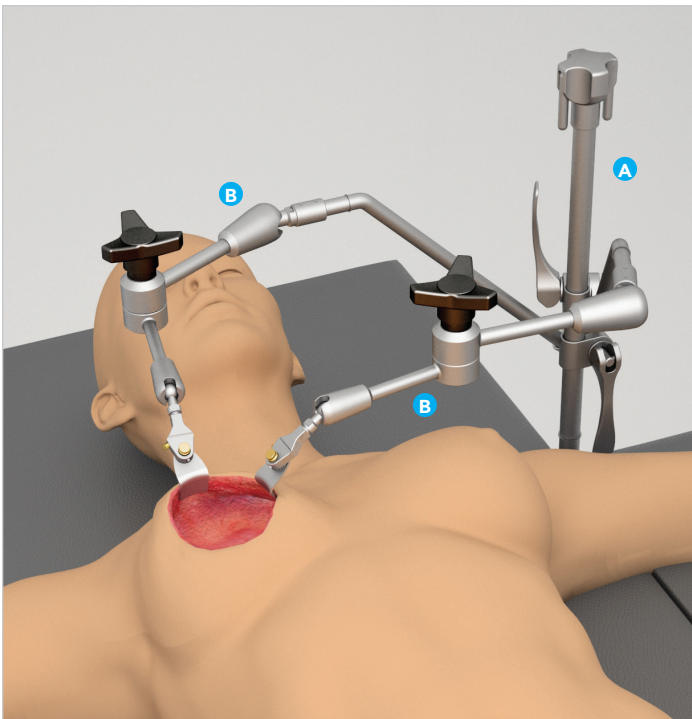

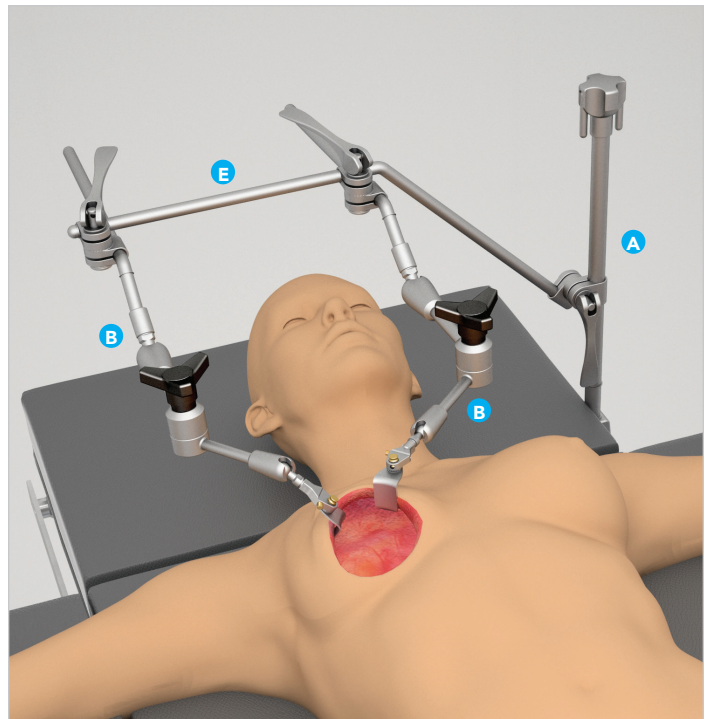


Mastectomy – Articulating Arm






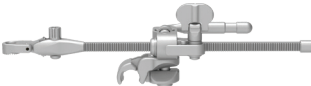


Articulating Arm (Alternate)




BREAST SYSTEM #SL82041

REF	IMAGE	PART #	ITEM DESCRIPTION	QTY	APPLICATION
Tray 1					
A		43903AC	Elite III Rail Clamp w/2 Cam Joints 18"	1	Attaches to the rail of the OR table and holds the articulating arm in position.
B		SL42191	S-Lock Articulating Arm Small 15 1/2"	2	Attaches to the rail clamp or crossbar and holds blades in place around incision.
		SL46126	Richardson 25mm x 32mm (1" x 1 1/4")	2	Muscle retraction for exposure with small or mini-incision.
		SL46422ET	Radio Richardson 50mm x 60mm (2" x 2 3/8")	2	Breast incision or abdominal wall retraction.
		SL46426ET	Radio Richardson 50mm x 100mm (2" x 4")	2	
		SL46432ET	Radio Richardson 50mm x 160mm (2" x 6 1/4")	2	
		SL46436ET	Radio Richardson 50mm x 200mm (2" x 7 7/8")	2	
		SL46436ET	Radio Richardson 50mm x 240mm (2" x 9 1/2")	2	
		SL46402	Cat's Paw 6mm x 38mm (1/4" x 1 1/2")	2	Skin retraction.
		SL46214	Murphy Sharp 4 Prong 25mm x 13mm (1" x 1/2")	1	Skin and shallow muscle retraction.
		SL46115	Malleable 13mm x 64mm (1/2" x 2 1/2")	1	Soft tissue retraction.
		SL46116	Malleable 13mm x 152mm (1/2" x 6")	1	Soft tissue retraction.
		SL46119	Malleable 25mm x 127mm (1" x 5")	1	Soft tissue retraction.
		50000G	Instrument Case (23" x 11" x 3 1/2")	1	Holds instruments for sterilization, storage, and transportation.

BREAST SYSTEM #SL82041 (Continued)

REF	IMAGE	PART #	ITEM DESCRIPTION	QTY	APPLICATION
Tray 2					
C		44615FH	14" Extension Arm with Fixed Hinge Joint	1	Attaches to the rail clamp and holds C-Frame in place.
D		44225	C Frame 13 1/2"	1	Alternate frame set up.
		SL42126WGLP	Cam II Clip-On Quick Angle Handle 10"	1	Attaches retractor blades to frame around the incision and angles 45 degrees up or down.
		SL45011CGP	Mini Micro-Adjustable II Clip-On Quick Angle 9"	1	
		60020	T-Handle	1	May be used on Quick Angle Handle to angle blades 45 degrees up or down.
		50000G	Instrument Case (23" x 11" x 3 1/2")	1	Holds instruments.
Tray 3					
		40020X	Lite Wand	1	Illuminates incision area.
		50000T	Instrument Case (11" x 7" x 3 1/2")	1	Holds light wand.

OPTIONAL ADD-ON CROSSBAR

REF	IMAGE	PART #	ITEM DESCRIPTION	QTY	APPLICATION
E		43900BC	24" Crossbar with 2 Cam III Joints (11" x 13")	1	Attaches to the rail clamp and holds the angled arms in place around the incision. May also support retractor handles.