


Knee Retractor Set

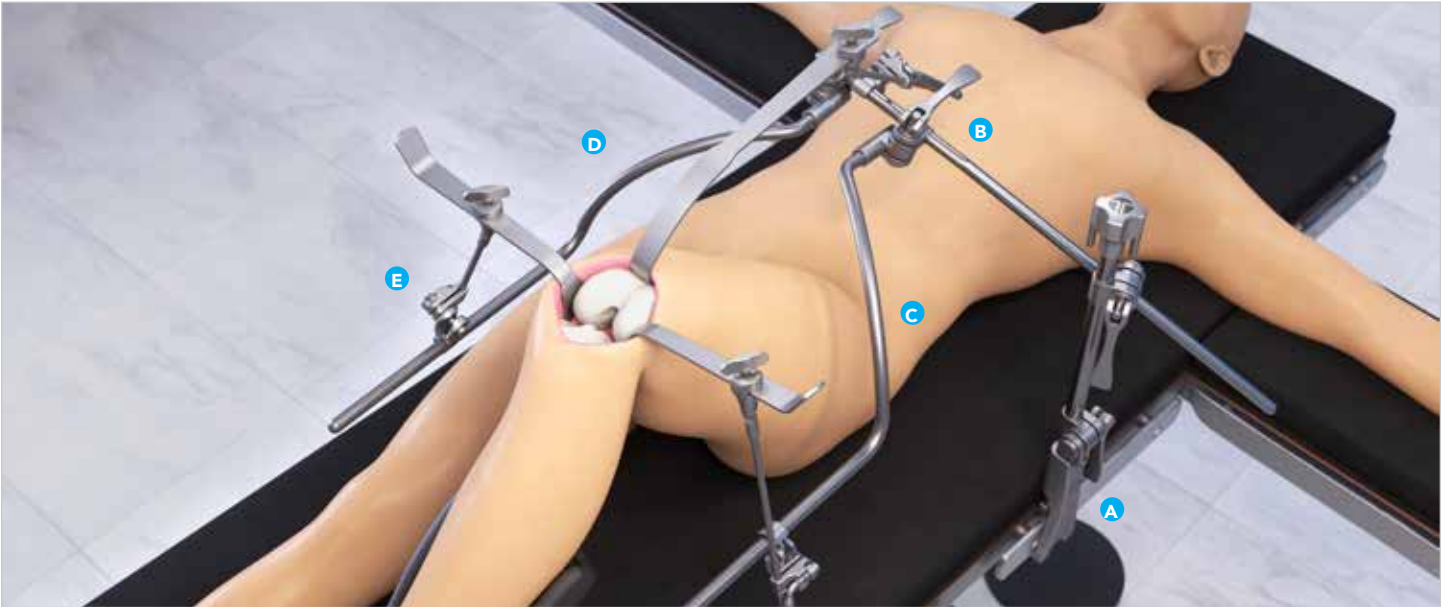

NOTE: Hospital's blades required for use.

HIP FRAME SMOOTH ARM KIT #82015

NOTE: Cleated arms are available for purchase if tether-only use is desired.

REF	IMAGE	PART #	ITEM DESCRIPTION	QTY	APPLICATION
Tray 1					
A		43902ACL	Elite III Rail Clamp with 1 Joint 16"	1	Attaches to the rail of the OR table and holds the crossbar in position.
B		43991	22" Crossbar with 2 Cam III Joints	1	Attaches to the rail clamp and holds the curved arms in place around the incision. May also support retractor handles.
		43646C	20" Upright Arm with Cam III Joint	1	Attaches to the rail clamp and may hold one curved arm in place..
C		44201	20" Curved Arm - 1	1	Attaches to the crossbar and holds the retractor handles and tethers in place around the incision.
D		44202	20" Curved Arm - 2	1	
		42103	Clip on Joint with Cleat	2	Optional attachment to Curved Arms if use of tethers is desired.
		42105	Hip Tethers 24" - 5 Pack (Latex Free)	1	Secures retractor blades to the frame.
E		42126BP	Cam II Clip-On Blade Adapter 8 1/2"	3	Secures retractor blades to the frame.
		42126BLP	Cam II Clip-On Blade Adapter 12"	1	
		50000G	Instrument Case (23" x 11" x 3 1/2")	1	Holds instruments for sterilization, storage, and transportation.
		41938	Rail Extender 20" with 1 Power Rail Clamp	1	Attaches to OR table rail to provide more attachment options for the rail clamp. Required to use Hip System with Hana® Table or leg positioning systems.