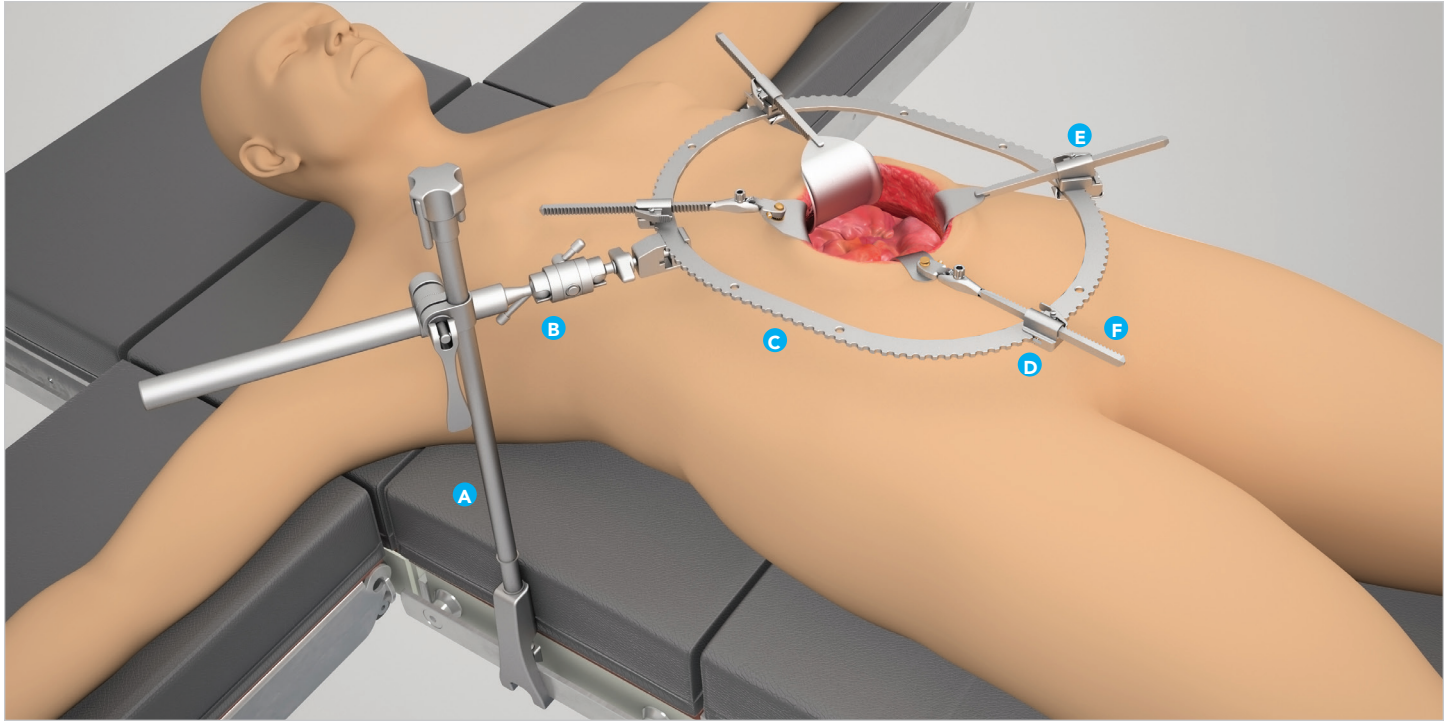








# RingTrack® General Abdominal


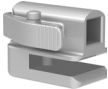







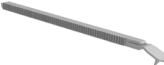



## RINGTRACK® GENERAL / VASCULAR / TRAUMA SYSTEM #SL82034

Fixed system also available.

REF	IMAGE	PART #	ITEM DESCRIPTION	QTY	APPLICATION
<b>Tray 1</b>					
A		47902ACL	RingTrack® Elite III Rail Clamp 16" with 1 Cam Joint	1	Attaches to the rail of the OR table and holds the horizontal bar in position.
B		47400	RingTrack® Articulating Horizontal Bar 18 1/2"	1	Attaches to the crossbar and holds the ring in place around the incision.
		47437	RingTrack® Horizontal Bar 14 1/4"	1	
		47418	RingTrack® Round Ring Medium Half 10 1/2"	2	Two are joined together to form a ring. Ring attaches to the horizontal bar and holds the retractors in place around the incision.
		47434	RingTrack® Straight Segment 8"	2	Two may be applied to two ring halves to create an oval ring
C		47412	RingTrack® Oval Ring Medium 10 1/2" x 14 1/2"	1	Attaches to the horizontal bar and holds the retractors in place around the incision.

**RINGTRACK® GENERAL / VASCULAR / TRAUMA SYSTEM #SL82034 (Continued)**

REF	IMAGE	PART #	ITEM DESCRIPTION	QTY	APPLICATION
		50000W	Instrument Case 24" x 16" x 3"	1	Holds instruments for sterilization, storage, and transportation.
<b>Tray 2</b>					
D		47440	RingTrack® Ratchet	2	Attaches to the ring and holds the retractor handle in place.
E		47442	RingTrack® Angling Ratchet	6	Attaches to the ring and holds the retractor handle in place. Allows the handle to be angled down.
F		SL45005G	RingTrack® Quick Angle Handle 8"	8	When used with the ratchet, attaches retractor blades to the ring and angles retractor blades as needed, 45° up or down.
		60020	T-Handle	2	May be used on Quick Angle Handle to angle blade 45° up or down.
		SL46160	Harrington 64mm x 152mm (2 1/2" x 6")	1	Organ or soft tissue retraction.
		SO46141C	Balfour [sm lips] 78mm x 73mm (3" x 3")	2	Costal margin or abdominal wall retraction.
		SL46154	Richardson 51mm x 76mm (2" x 3")	2	Abdominal wall retraction.
		SL46150	Richardson 51mm x 102mm (2" x 4")	1	
		SL46151	Richardson 51mm x 127mm (2" x 5")	1	
		SL46152	Richardson 51mm x 152mm (2" x 6")	1	
		SL46112	Malleable 51mm x 165mm (2" x 6 1/2")	2	Soft tissue retraction.
		SL46120B	Malleable 38mm x 152mm (1 1/2" x 6")	1	
		SO46124A	Malleable 76mm x 254mm (3" x 10")	1	
		47115W	RingTrack® Gelpi 1/2"	1	
		50000G	Instrument Case (23" x 11" x 3 1/2")	1	Holds instruments for sterilization, storage, and transportation.